

The Eliot Patient Chair Range



The Eliot Patient Chair Range

RMA are pleased to re-launch the Eliot chair formerly manufactured by James Spencer.

Using our state-of-the-art production facility based in South Wales, RMA set out to improve the existing design through the development of single action leg height adjustment which eases the process of adjusting the chair height for the caregiver. The transit feature from the original chair has been retained and improved by reducing the complexity of the bolt fixings whilst still retaining the ability to service the castors if and when required.

Fixed Arm

Code 5770



Flip-Arm Standard

Code 5771



Flip-Arm Bariatric

Code 5775



The height adjustment caters for pressure area care on the heels and sacral area in varying patient heights. The transit castor enables easy transportation of the chair.




The open architecture of the chair allows safe and easy patient evacuation. The gap between seat base and back enables sling access for hoisting.




Colour Options



Upholstery Specifications

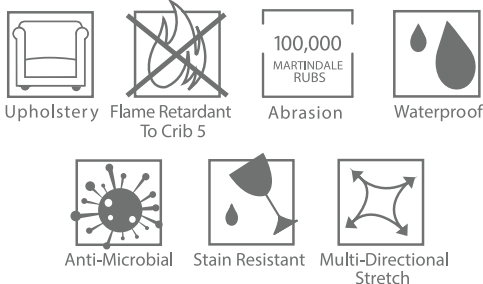

Head Position, Neck and Spine Alignment

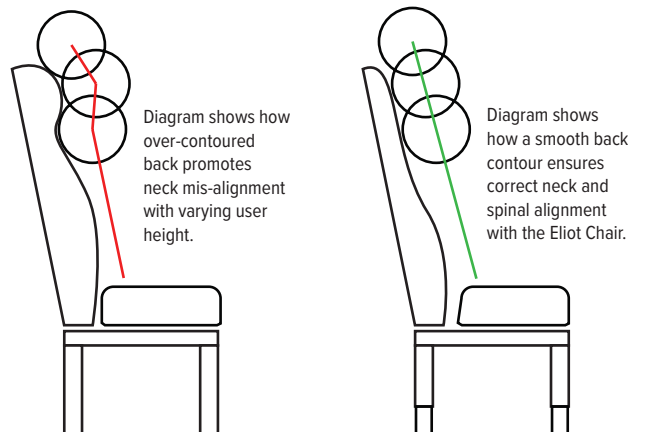

The Eliot Bed Side Chair

Code 5770

Specifications

- ▶ Powder Coated Steel Frame
- ▶ Lacquered Wood Arm Pads
- ▶ Height Adjustable
- ▶ Flame Retardant Upholstery
- ▶ Seat Pad Gap for hoist sling access
- ▶ Transit Castors
- ▶ Pressure Relieving Memory Foam
- ▶ Floor to Seat Pad Height: 400 to 500mm
- ▶ Overall, Height: 1140 to 1240mm
- ▶ Overall Depth: 720mm
- ▶ Overall Width: 600mm
- ▶ Seat Pad Depth: 460mm (490mm to padded back)
- ▶ Seat Pad Width: 460mm (470mm between arms)
- ▶ Weight of product: 21Kg
- ▶ Maximum user weight: 190Kg (30St)
- ▶ The chair meets BS 7176:2007
- ▶ ISO 9001:2015 Accredited



Smooth back contour ensures correct neck and spinal alignment.

The Eliot Flip Arm

Code 5771

Specifications

- ▶ Powder Coated Steel Frame
- ▶ Flip-up Arms
- ▶ Height Adjustable
- ▶ Flame Retardant Upholstery
- ▶ Seat Pad Gap for hoist sling access
- ▶ Transit Castors
- ▶ Pressure Relieving Memory Foam
- ▶ Floor to Seat Pad Height: 400 to 500mm
- ▶ Overall, Height: 1160 to 1260mm
- ▶ Overall Depth: 720mm
- ▶ Overall Width: 600mm
- ▶ Seat Pad Depth: 460mm (490mm to padded back)
- ▶ Seat Pad Width: 460mm (470mm between arms)
- ▶ Weight of product: 21Kg
- ▶ Maximum user weight: 190Kg (30St)
- ▶ The chair meets BS 7176:2007
- ▶ ISO 9001:2015 Accredited



Height Adjustable Legs



Flip-Up Arms



The Eliot Bariatric

Code 5775

Specifications

- ▶ Powder Coated Steel Frame
- ▶ Flip-up Arms
- ▶ Height Adjustable
- ▶ Flame Retardant Upholstery
- ▶ Seat Pad Gap for hoist sling access
- ▶ Transit Castors
- ▶ Pressure Relieving Memory Foam
- ▶ Floor to Seat Pad Height: 400 to 500mm
- ▶ Overall, Height: 1140 to 1280mm
- ▶ Overall Depth: 760mm
- ▶ Overall Width: 680mm
- ▶ Seat Pad Depth: 455mm (550mm to padded back)
- ▶ Seat Pad Width: 680mm
- ▶ Weight of product: 30Kg
- ▶ Maximum user weight: 381Kg (60St)
- ▶ The chair meets BS 7176:2007
- ▶ ISO 9001:2015 Accredited

Flip-Up Arms



← - - - Extra Wide Bariatric Width - - - →
680mm



Seat can be adjusted forward increasing the distance between the seat and the seat back. This customisation accommodates the variance of seat depth required.





York Road
Bridgend Industrial Estate
Bridgend
CF31 3TB
Tel. 01656 674488
Fax. 01656 674499
Email. sales@romamedical.co.uk
Web. www.romamedical.co.uk

